



Inc WIN device settings - Version 1.4.1 (2023)

Window Incubator / WIN

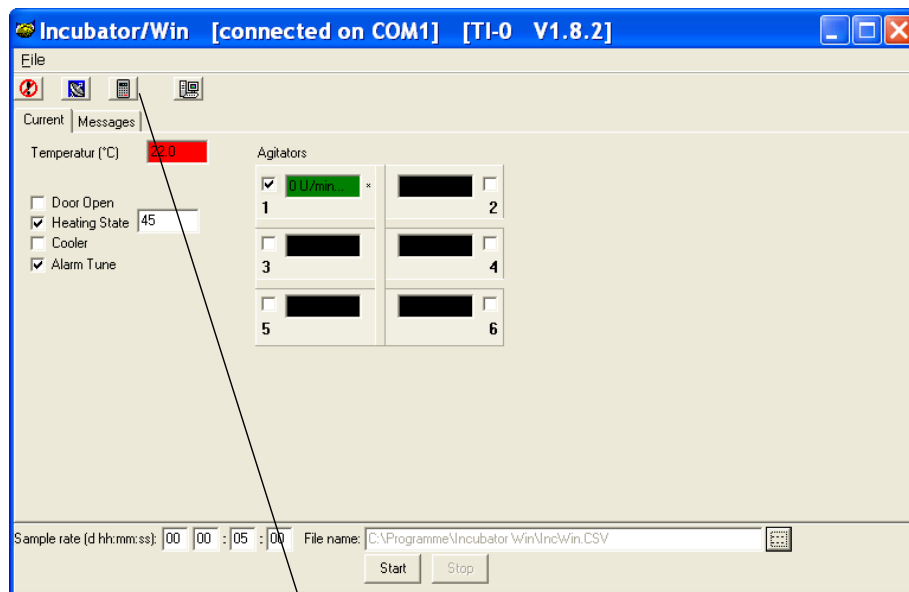

Incubator Configuration Basic setting

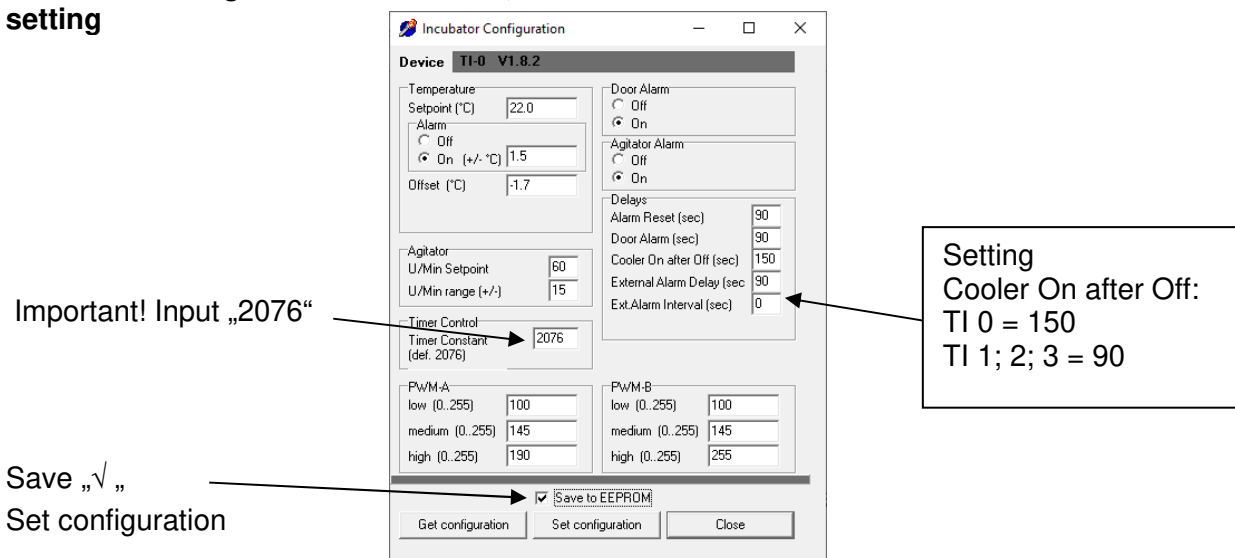

PWM A+B settings depending on device type:

		TI 0	TI 1	TI 2	TI 3
Fan speed low	PWM-A low	100	100	100	120
Fan speed middle	PWM-A medium	145	145	145	145
Fan speed high	PWM-A high	190	190	175	200
Heating low	PWM-B low	100	100	100	120
Heating middle	PWM-B medium	145	145	145	145
Heating high	PWM-B high	255	255	255	255